


PRO-TEC® III

SYRINGE SHIELD

The Pro-Tec® III Syringe Shield reduces hand exposure and is convenient to use. The barrel of the shield is constructed of .08" (2 mm) thick tungsten that will reduce radiation exposure from Tc-99m by 94%. A 5.05 density (1.5 mm lead equivalency) lead glass window provides protection and visibility. The interior of the barrel behind the glass has a white coating that improves the view of the syringe and contents. A bevel around the lead glass helps protect it from scratching or breaking.

The Safe-T-Lock design provides increased protection by allowing minimal handling. The Safe-T-Lock design grips and secures the syringe upon insertion and will release the syringe at the press of a button. Disposing of used syringes is easy; invert the syringe shield over a sharps container, press the release button and the syringe freely disengages.

Pro-Tec® III Syringe Shields accommodate the standard sized 1 cc, 3 cc, 5 cc, 10 cc and 20 cc syringes.

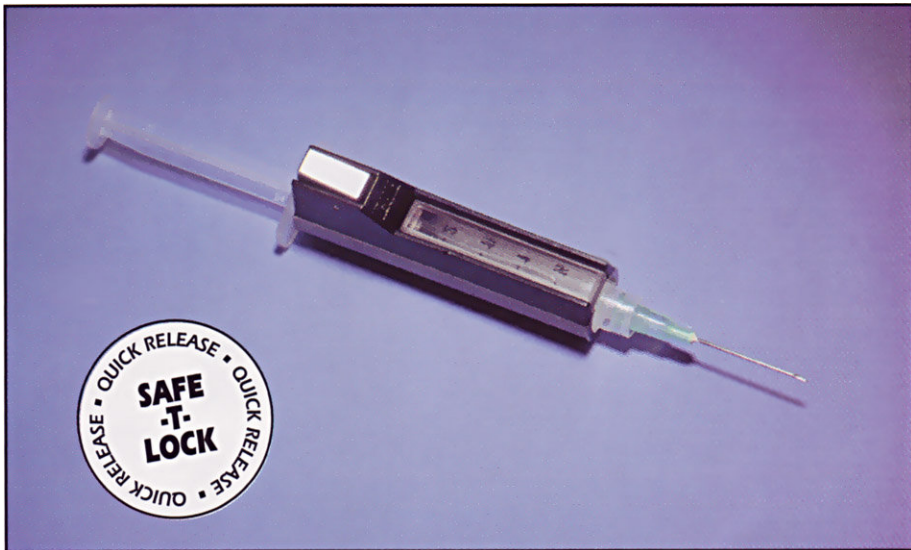

Longer design accommodates most syringes

FEATURES:

- *Tungsten shielding*
- *Lead glass window*
- *Lightweight*
- *Fits most disposable syringes*
- *Safe-T-Lock design reduces exposure with faster handling*

Pro-Tec® III Syringe Shields:

007-734	Syringe Shield, 1 cc	\$145.00
007-735	Syringe Shield, 3 cc	145.00
007-736	Syringe Shield, 5 cc	160.00
007-738	Syringe Shield, 10 cc	185.00
007-747	Syringe Shield, 20 cc	395.00

Pro-Tec® III Syringe Shields (for international syringes):

007-755	Syringe Shield, 2.5 cc	\$145.00
	<i>For Terumo syringe</i>	
007-745	Syringe Shield, 3 cc	145.00
	<i>For BD Plastipak and Terumo 3mil syringes</i>	
007-746	Syringe Shield, 5 cc	160.00
	<i>For Once syringes</i>	

Prices subject to change without notice.

www.jgravengaard.com 2007-12-04



Safe-T-Lock design quickly engages or releases syringe.